



DIMAGRELLE X·TONE

**SLIMS
MODELS
FIRMS**

BIOSTIMULATING IMPULSE
EQUIPMENT FOR BODY
REMODELING, TONING
AND MUSCLE FIRMING.



DIMAGRELLE X-TONE

All the actions normally carried out with physical activity can now be performed at the institute with the new Dimagrelle X-Tone stimulation device.

Slimming, cellulite reduction, increase in tone and firming of the tissues of the whole body in a practical, fast and effortless way.

Dimagrelle X-Tone carries out stimulation treatments suitable for every morphology: Body Shaping & Muscle toning and Localized Firming & Lipomodeling.



BODY SHAPING & MUSCLE TONING



QUICKLY SLIMMING, DRAINING AND ACTING ON CELLULITE


Total Body stimulation treatments aim at general metabolic reactivation, body remodeling, reduction of cellulite blemishes and vascularization, drainage and tissue toning.



DIMAGRELLE TOTAL BODY TREATMENTS

In this type of treatment, a special muscle stimulation with variable impulses is delivered, which activates the muscle pumps and obtains drainage and toning of the tissues: this specific stimulation allows the simultaneous delivery of the Dimagrelle products and enhances their absorption and specific properties. The treatment is performed in a special electric blanket that enhances both the bio-stimulating action and the absorption the products. A complete and very satisfying treatment!

LOCALIZED FIRMING & LIPOMODELING

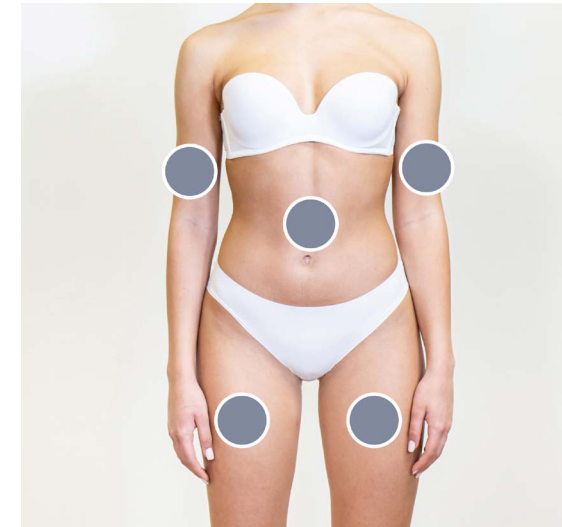

FIRMING AND SHAPING THE CRITICAL AREAS

Localized stimulation treatments aim to activate and firm the muscle groups of critical areas such as abdomen, arms, buttocks and thighs.



PROGRAMS FOR EVERY NEED

Depending on the customer's needs and on the area to be treated, the working program can be customized by choosing between 5 different types of waves. These activate the metabolism and muscle work, necessary to obtain localized remodeling and local lipolytic-slimming and firming effects of the various body areas in a few sessions.

BODY SHAPING & MUSCLE TONING

Total Body combined stimulation treatment

This type of treatment is suitable for acting on water retention, cellulite blemishes and tissue atony. It produces a marked effect of reduction and remodeling of the shape, visible from the first session.



The treatment combines 2 types of electrical stimulation that perform a synergistic action:

- Delivery of active substances contained in specific products.
- Stimulation of muscle pumps to activate the vascular system and tone the tissues.

To amplify the slimming and toning effects, the treatment is performed with the combination of an electric blanket with 3 independent zones of heating: torso / abdomen / legs. The heating generated, in combination with the products and electrical bio-stimulations, allows for effective drainage, acceleration of metabolisms, reduction and perfect remodeling of the shapes.

Total Body programs allow the specific treatment of the whole body or of more concentrated areas. And depending on the needs of the treated subject, they can be performed, even alternately, with localized muscle stimulation treatments.

The result will progressively lead to vascular reactivation, the decrease of localized accumulations and the redefinition of the shapes, with an increase in the general tone of the tissues.

A COMPLETE REMODELING TREATMENT



BIOSTIMULATION



PLANTS
EXTRACTS



VEHICULATING
STIMULATION




THERMO-SLIMMING
ACTION

BODY SHAPING & MUSCLE TONING

The combined action of heat, medicinal herb extracts and bio-stimulating impulses allows for a mix of actions that remodel, tone, drain and activate the tissues simultaneously, while providing a pleasant and relaxing treatment lasting about an hour.



- 1** Apply the **Dren concentrated gel** or **Snell concentrated gel** on the various areas of the body, as indicated in the treatment sheets.

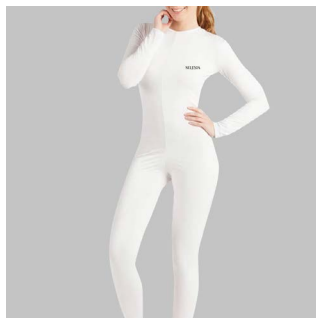

- 2** Put on the **elastic suit** previously soaked in the Dimagrelle Osmodynamic solution with an osmotic-purifying and draining action.



- 3** Apply the **conductive bands** according to the diagrams indicated.



- 4** Connect the **electrodes to the conductive bands** respecting the indicated polarity.



- 5** Wrap the client in the **electric blanket**, differentiating the temperature in the various parts of the body, as needed, and start the biostimulation treatment.



- 6** End the **treatment** with the application of the Intensive Cellulite Emulsion with a short final massage.



LOCALIZED FIRMING & LIPOMODELING

Firming Biostimulating Treatment

Localized programs for muscle toning and firming allow the treatment of laxity and fatty deposits in critical areas. The different type of waves allow to customize the work for every need and to find and maintain a firm and harmonious body and a slender silhouette. The muscle contractions produced firm and predispose the body to greater energy consumption and, therefore, to a decrease in localized accumulations. This type of treatment can be completed and made even more performing with the combination of an electric blanket and Body Concept toning products.



MUSCLE
STIMULATION



FIRMING AND
INCREASING STRENGTH



REDUCTION OF ADIPOSE
ACCUMULATIONS

LOCALIZED FIRMING & LIPOMODELING

The critical points are treated with specific stimulations that activate the deep muscle groups, obtaining an effective lipolytic and firming action that redefines the shapes. Each localized treatment lasts from 30 to 60 minutes per area.



It firms the muscles (biceps and humeral triceps) and plumps the tissue reducing laxity.



It firms, sculpts and reduces localized fat. Specific action on the muscles of the rectus abdominis and mm.



It firms and reshapes the buttocks, giving a firm and slender appearance.



It firms and remodels the quadriceps with lipolytic and anti-cellulite action.



It firms and tones the ischio-peroneo-tibial muscles giving a firm and slender effect.

LOCALIZED FIRMING & LIPOMODELING


60-minute treatment that combines biostimulation, heat and specific toning products: a perfect blend to obtain an effective metabolic action. This mix of combined actions gives pleasure, completeness and value to the basic treatment.




1 Prepare and turn on the electric blanket opened on the bed and cover it with a protective paper sheet.



2 Apply the electrodes for muscle groups after having suitably sprinkled them with the Biostimulation Gel.

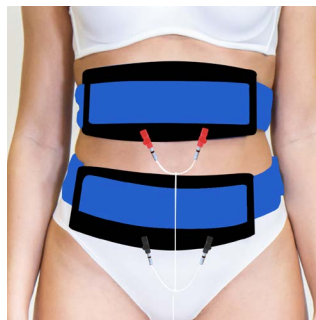


3 Connect the cables and plugs to the electrodes and close the cartene sheet and the electric blanket.

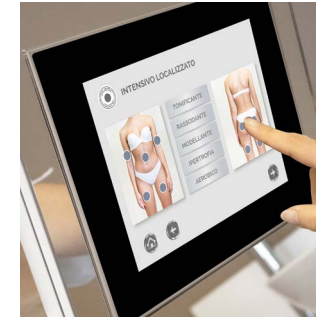

4 Choose the stimulation program and start the treatment. Duration of programs 30 '.



5 Open the Electric Blanket and remove the electrodes.

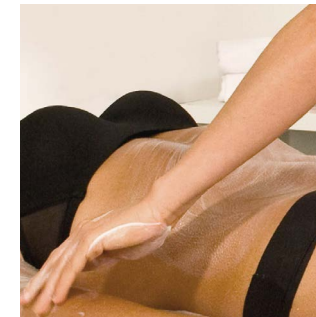

6 Apply the combined products and close the electric blanket, extending the session for about 25 minutes.



7 End the treatment by allowing any residues of the product already applied to penetrate with a short massage.

Programs for every need

Five biostimulation programs allow to carry out personalized treatments for every need.

- 1. Toning
- 2. Firming
- 3. Modeling
- 4. Hypertrophy
- 5. Aerobic

1 TONING

Program suitable for the treatment of tissues with water retention, slight atony and localized cellulite. Suitable for the first toning sessions, to prepare the muscles before the most intense firming work.

2 FIRMING

To be used as the main muscle firming program.

3 MODELING

To be used as a "sculpting" program of the muscles already treated with the firming program.

4 HYPERTROPHY

It increases the diameter and capacity of muscle fibers by restoring muscle volume and strength.


5 AEROBIC

It increases the muscle tone and improves circulation and oxygenation of the treated muscles. It helps to delay fatigue and loss of tone.



DISPLAY

A convenient display allows you to select the most suitable type of program for each imperfection and area, with indication of the points of application of the electrodes in the various areas, in a practical and intuitive way.

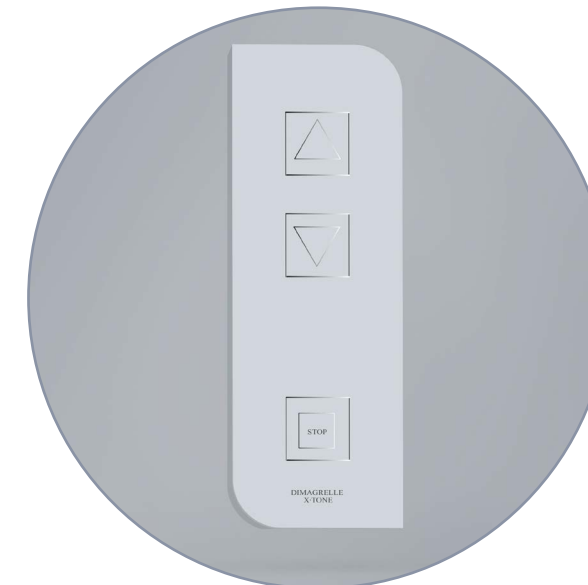

REMOTE CONTROLLER

A practical remote control allows the client who undergoes the treatment to adjust the intensity of the stimulation in an optimal way and in an absolutely autonomous way, thus combining effectiveness and pleasure at the highest levels.




DIMAGRELLE X-TONE

A complete technology that brings numerous advantages to the beauty center!

- **TOTAL BODY AND LOCALIZED TREATMENTS**
Dimagrelle X-Tone is suitable for any body imperfection and is used effectively on the whole body and on the most critical localized parts such as arms, abdomen, back, thighs and buttocks.
- **FAST AND EFFECTIVE RESULTS**
Ideal for those who want quick and persistent results because it reactivates skin metabolisms, reduces cellulite and tones relaxed muscles.

- **PLEASANT AND RELAXING TREATMENT**
This non-invasive Bio-electrostimulation treatment is well tolerated by all people and the client himself can adjust the intensity of the stimulation.
- **HYGIENE GUARANTEE**
Perfectly hygienic, thanks to the practical washable and reusable elastic bodysuit and the use of disinfectable electrodes: maximum efficiency in compliance with anti-Covid 19 regulations.
- **WORK WITHOUT AN OPERATOR**
It allows you to perform modeling treatments without the constant presence of the operator and therefore allows you to have a higher gain.
- **PROFITABLE**
With a low cost for the beautician in terms of products, consumption and personnel, Dimagrelle X-Tone allows to offer effective and complete high-profitability treatments and programs.

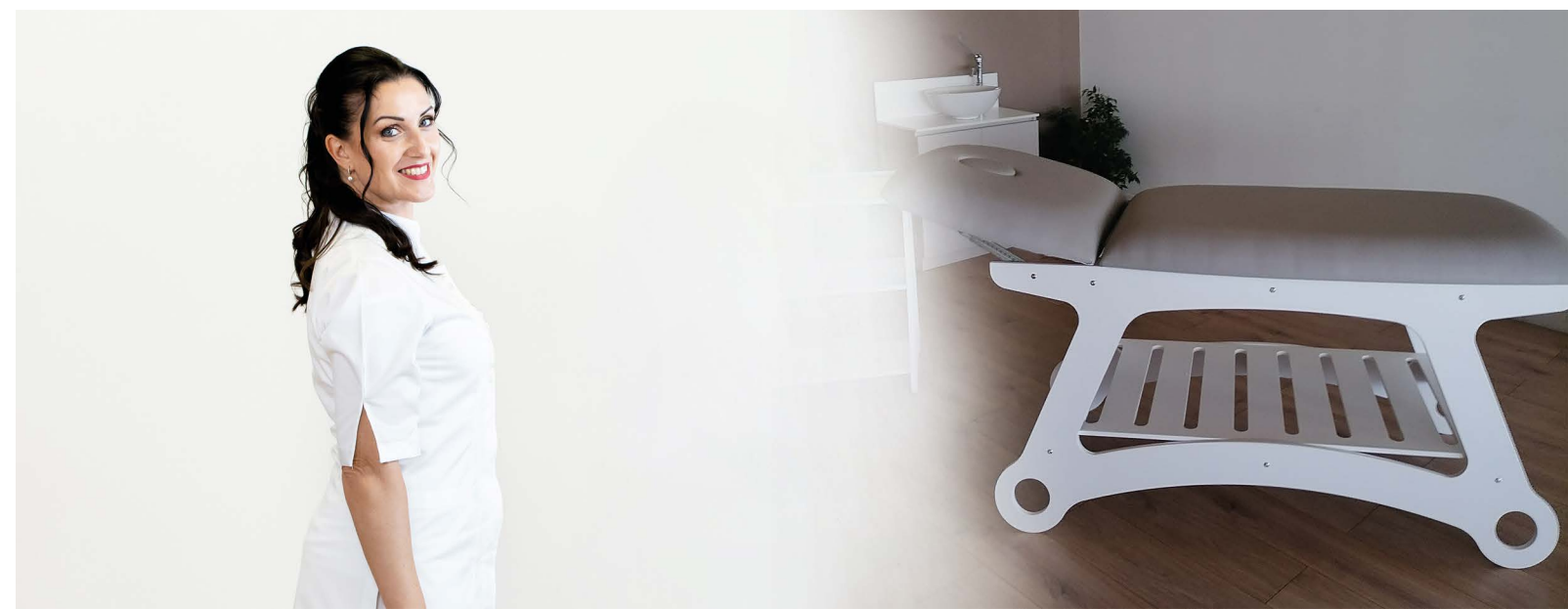



DIMAGRELLE
X·TONE



SELÈNIA
ITALIA

Via Vicarese 81 - 56012 Calcinaia (PI) - Tel. +39 0587 48 99 77 - selenia@seleniaitalia.it
WWW.SELENIAITALIA.IT